User Manual

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Document No

Doküman No

TF03-IFU03-06

Issue Date

03.06.2022

Revision No

01

Revision Date

11.08.2022

**SIMCOE AND ANTERIOR CHAMBER STABILIZATION
CANNULA**
**SIMCOE ve ANTERIOR CHAMBER STABİLİZASYON
KANÜLÜ**

Company Name

Şirket Adı

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Web

: www.merantip.com

Product Structure and General Features:

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The Simcoe and Anterior Chamber Stabilization Cannula is designed for the use in the cataract surgeries to maintain anterior chamber. Cannula tip is made of 304 Stainless Steel. The connection area at the back of the cannula is made of polycarbonate raw material. The products are available for single use and sterile. The cannulas to be used may differ depending on the operation area and the surgeon's choice. Models are available in various diameters and sizes in order to be suitable for every eye structure.

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Product Models:

Ürün Modeli:

Model and size information of the Simcoe and Anterior Chamber Stabilization Cannula products manufactured by us are shown in **TF03.TB.01 Annex-1 Product Models Table** with pictures.

Annex-1 Product Models Table with pictures.

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Revision Date
11.08.2022

Intended Use:

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Simcoe Cannulas are disposable and usually used in cataract surgeries. Anterior Chamber Stabilization Cannulas are using especially by Anterior Chamber Stabilization operations with phacoemulsification devices.

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Indications:

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Simcoe Cannulas are used to provide a method for irrigation and aspiration without the need for a phaco machine. Anterior Chamber Stabilization Cannulas provide irrigation for chamber maintenance during anterior surgeries.

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Contraindications:

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It should not be used other than ophthalmic surgery.

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It should not be used on infected eyes.

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Complications:

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If the tip of the Cannula is not sharp, it may damage the tissue of the eye.

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After operation mild chemosis and redness can occur.

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Intraoperative serous occurrence.

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Side Effects:

Yan Etkiler:

There are no known side effects of the product.

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Product Users:

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The product should be used by specialist physicians who have completed the necessary medical training. Users must have operational experience.

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Instruction for Use:

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Simcoe and Anterior Chamber are restricted for use by surgeons trained in ophthalmic surgery.

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➤ The cannulas should be opened from the marked are without disruption of the sterility.

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➤ Cannula should be connected to the Equipment by a sterile person according to usage of the product.

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➤ The cannulas are applied in anterior chamber of eye for irrigation and aspiration purposes.

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➤ Dispose of used products in accordance with established hospital protocol for hazardous waste.

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Factors That Can Affect the Success of The Operation:

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Protect the product against direct sunlight during shipment and storage.

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Keep the product dry in any environment.

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Store the products between 5 ° C and 35 ° C temperature and maximum 70% RH humidity.

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Medical Equipment Used with Products:

Medical Equipment Used with Products:

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Revision Date
11.08.2022

Packaging, Labeling and Handling:

Packaging material are manufactured according to EN ISO 11607-1 and EN 868-6 standard, tear, puncture and high resistance to liquids, which maintains sterility is the long packaging material. It sealed with heat cut in desired dimensions. Leak-proof barrier allows the product to remain sterile.

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
Revision Date

11.08.2022

The packaging sample of our products is as seen in the photos below.

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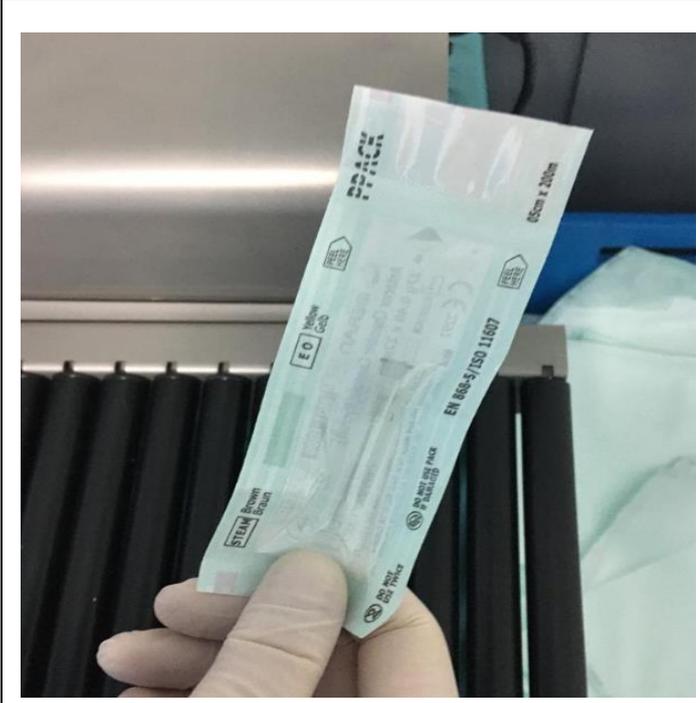Revision Date

11.08.2022

Symbols Prepared According to EN 15223-1 Standard and Their Meanings

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